

PFT™ Powder Flow Tester

...affordable testing for powder characterization

The PFT Powder Flow Tester brings quick and easy analysis of powder flow behavior in industrial processing equipment. Evaluate powder discharge from storage containers. Use as QC check for incoming materials. Rapidly characterize new formulations for flowability and adjust composition to match flow behavior of established products.

Choice of Test Options:

- Flow Function
- Time Consolidated Test with Flow Function
- Wall Friction
- Bulk Density

Choice of Flow Function Tests:

- Standard 5 Point (16 minutes)
- 2 Point (10 minutes)
- Time Consolidation

Real Time Clock Displays:

- Test Step
- Remaining Time to Completion

Shearing Algorithm Captures:

- Peak Stress Value
- Subsequent Stable Stress Value

Data Output Flow Function:

- Flow Index for Powder Flowability
- Arching Dimension
- Rat-hole Diameter
- Hopper Half Angle
- Gravity Chute Angle
- Wall Friction Angle
- Bulk Density

Compact design with small footprint
Tester fits conveniently on workbench
Depth: 15inches / 38cm
Width: 14inches / 36cm
Height: 27inches / 69cm



Developed in association with
The Wolfson Centre
for Bulk Solids Handling Technology
at the University of Greenwich, England.

What's Included?

Instrument
Powder Flow Pro Software
with USB Cable

What else is needed?

Choose one or both:

- Standard Volume Accessory Kit
230cc Trough & 33cc Vane Lid
- Small Volume Accessory Kit
38cc Trough & 5cc Vane Lid

Each of the above include:

- Wall Friction Lid
304 s/s simulated 2B finish
- Outer Catch Tray
- Inner Catch Tray w/ Scraper Tool
- Powder Scoop
- Cleaning Brush

Optional Accessories

- Wall Friction Lids in Mild Steel
22-28RA, Tivar 88 or special order
- Temperature Probe
- Humidity Sensor
- Sieve Kit Standard or Small Volume
- Powder Flow Demonstration Tool




Small Volume Vane Lid
.795-13.252 kPa

Standard Volume Vane Lid
.289-4.819 kPa

Powder Flow Pro Software Included


Operation and control of the Powder Flow Tester is accomplished with Powder Flow Pro Software.

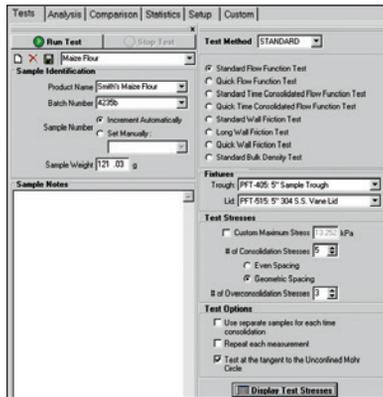


Main screen provides choice of basic tests:

- Flow Function
- Wall Friction
- Time Consolidated Flow Function
- Bulk Density

Custom test options available

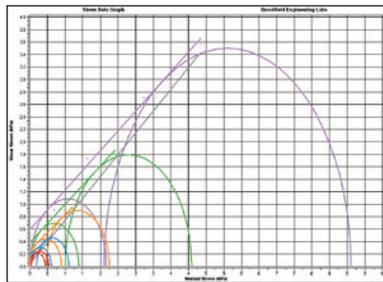

Stress data output screen captures "normal stress" and "shear stress" values and plots data in graphical format (calculates Mohr Circle Failure Loci).

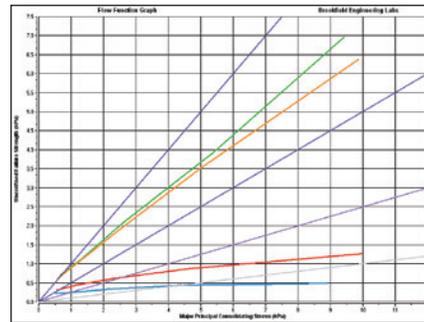


Flow Function test produces graphs of powder flow behavior which show:
 Unconfined Failure Strength vs. Major Principal Consolidating Stress
 Arching Dimension vs. Principal Consolidation Stress
 Rat-hole Diameter vs. Powder Fill Level
 Bulk Density (Fill and Final)



Close-up View of Vane Lid used for Flow Function Test.



Close-up View of Wall Friction Lid for wall friction and bulk density test

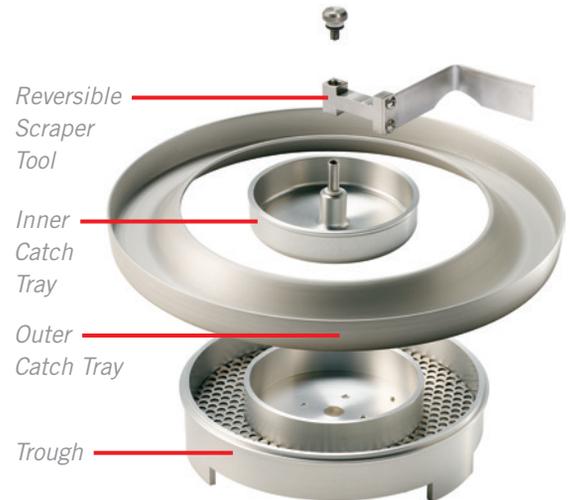

Outer and Inner Catch Trays with Scraper Tool for Sample Preparation in Trough

PFT POWDER FLOW TESTER SPECIFICATIONS	
Load for Vertical Axis Compression:	7 kg — Accuracy $\pm 0.6\%$ FSR
Axial Speeds:	0.1mm/second up to 5mm/second
Distance:	Accuracy $\pm 0.3\text{mm}$
Torque:	$\pm 7.0 \text{ N}\cdot\text{m}$ — Accuracy $\pm 1.2\%$ FSR
Trough Rotational Speeds:	1 revolution/hour (RPH) up to 5 RPH
Temperature Sensing:	-20°C to 120°C^*
Humidity Sensing:	10% to 95% RH $\pm 5\%^\dagger$
Dimensions (wxdxh):	(cm) 36.2 x 39.7 x 67.6 (in) 14 $\frac{1}{4}$ x 15 $\frac{5}{8}$ x 26 $\frac{5}{8}$
Weight:	34 kg (75 lb)

* Requires Part No. DVP-94Y

† Requires Part No. PFT-607Y

MINIMUM COMPUTER SPECIFICATIONS FOR POWDER FLOW PRO SOFTWARE	
2GHz processor with 1 GB of RAM and 30 MB hard drive space available	
1024x768 video resolution with 128 MB of graphics memory	
Windows 7, 8 or 10 (32 and 64 bit) with one USB or RS-232 port	