# CT3 Accessories

FOR A WIDE RANGE OF TESTS.

While many CT3 accessories have multiple applications, some are particularly useful for specific industries. The following color coded icons are used here to identify these industries.

- FOODS
- COSMETICS
- PHARMACEUTICALS
- MATERIALS & PACKAGING
- **D** DEVICES MECHANICAL



TA-KF @

Kieffer Dough and Gluten Extensibility Fixture quantifies maximum force and distance needed to break sample. Fixture Base Table required.




TA-JPA 📵 🐠

Junior Punch Fixture is for punching through flat samples; 12.7mm max. diameter probe. Hole in fixture is 14mm. Rotary Base Table required.



TA-TPB @

Three Point Bend Fixture is used with TA7 blade from general probe kit. Fixture Base Table required. Small scale version TA-JTPB is available.

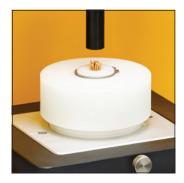

TA-DSJ **3** 

Dough Stickiness Fixture is standard test for measuring dough stickiness; important for processing raw dough. Fixture Base Table required.



TA-DE @

Dough Extensibility Fixture for holding sheet of raw dough or flat bread to measure breaking point of stretched sample. Fixture Base Table required.



TA-CTP **(3**)

Compression Top Plate for applying uniform compression forces on samples up to 4x6 inches (10x15cm) Fixture Base Table required.



TA-VBJ 🙃

Volodkevich Bite Jaws for testing bite force of meat products using shear cutting-test. Fixture Base Table required.



TA-PTF **6 0** 

Pizza Tensile Fixture quantifies cooked pizza firmness by measuring the tensile force and deformation distance to break sample.



TA-FMBRA **3** 

Standard dough pot set for preparing dough samples and measuring dough firmness.



TA-AACC36 @

AACC spec probe for measuring bread firmness and performing texture profile analysis (TPA). Fixture Base Table required.



TA-SBA-WB-1 G

1mm Shear Blades for cuttingshear test: meat, fish, sausage, etc. TA-SBA-WB-3 option for 3mm blades.

Fixture Base Table required.



TA-SFF @

Spaghetti Flexure Fixture quantifies flexure characteristics of uncooked spaghetti and other dry pastas.



#### TA-OC-002 **3**

Ottawa Cell (447cc) for bulk compression to determine hardness and crispness of cereals, vegetables or fruits. Fixture Base Table required.



### TA-PFS **3**

Pasta Firmness and Stickiness Fixture measures the firmness and stickiness of uncooked pasta. Fixture Base Table required.



### TA-PFS-C **3**

Pasta Firmness and Stickiness Fixture measures shear strength when biting pasta and like products. Fixture Base Table required.



### TA-KSC-002 **3**

Kramer Shear Cell - Sharp to measure shear force of small composite samples such as grapes, corn and beans. Fixture Base Table required.



### TA-CSF @

Circular Support Fixture provides support for round samples and retains any potential fluid expressed during the test. Fixture Base Table required.




# TA-WSP 6

Wire Shear Plate cuts through the sample. Good for products with significant stickiness like cheese and butter. Fixture Base Table required.



### TA-CEF **(3**

Cheese Extensibility Fixture measures the extensibility of molten cheese sample to breaking point.

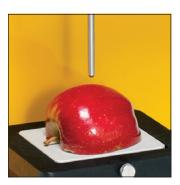

#### TA-MTP **()**

Magness-Taylor Probes for puncture test to measure hardness of fresh fruit and vegetables.

Fixture Base Table required.



### TA10 **3**

GMIA & GME probe and spec Bloom bottle TA-GBB-2 sold in package of twelve bottles.



### TA-MCF **G**

The Multiple Chip Fixture is used for testing the penetration or firmness of multiple chips / french fries.

Fixture Base Table required.



# TA-CJ 🧿

Confectionary Fixture for holding candies and similar products for penetration testing.


Fixture Base Table required. Probe not included.

# **CT3** Accessories

FOR A WIDE RANGE OF TESTS.

While many CT3 accessories have multiple applications, some are particularly useful for specific industries. The following color coded icons are used here to identify these industries.

- FOODS
- COSMETICS
- PHARMACEUTICALS
- MATERIALS & PACKAGING
- **D** DEVICES MECHANICAL

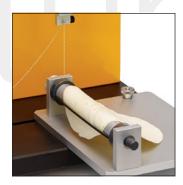

#### TA-TRF **(3**)

Tortilla Roll Up Fixture evaluates changes in corn tortilla texture per AACC technical paper by measuring the force to roll up a tortilla.



TA-CKA **3** 

Craft Knife Adapter cuts cleanly into and through material with minimum deformation of the sample.



TA-52 MOHRS **G** 

Shear Blade used for cutting tests, especially meat, poultry, fish or similar products. Fixture or rotary base table recommended.



TA-MP **(3** (0)

Mesh Probe quantifies the consistency of products such as mayonnaise and yogurt.



TA-BEC 📵 🐠

Back Extrusion Cell for measuring consistency of applesauce, pudding, yogurt and similar products. Rotary Base Table required.



TA-AVJ 🙃 🐠

Adjustable Vice Fixture for holding small samples for puncture test. Good for jelly beans, gum drops, etc. Rotary Base Table required.



TA-JMPA 🙃 😉

Multiple Probe Assembly consisting of nine 3mm probes and base plate designed to hold nine small samples of irregular geometry. Base Table required.



TA-DEC 🙃 😉 😉

Dual Extrusion Cell for either forward or back extrusion of fruit puree, pudding, yogurt or similar products. Fixture Base Table required.



TA-STF @ P G

Spread Test Fixture quantifies the spread force of a material. Comes with 1 male cone probe, 5 five samples cups and 1 sample cup holder. Base Table required.



TA-HCF (9)

Hair Combability Fixture measures the effect of hair dye, shampoos and conditioners on the combability of hair.



TA-EP 🔞 🐠

Eve Pencil Test Fixture measures hardness of cosmetic pencil tips for eye- or lip-lining products and can also be used for artistic type pencil tips. Fixture Base Table required.



TA-LC <sup>(3)</sup>

Lipstick Cantilever Test Fixture allows imitative tests on lipstick and similar products to quantify strength of product. Fixture Base Table required.



TA-TE @ @

Tube Extrusion Jig measures the force needed to squeeze cream, paste or ointment out of a tube. Base Table required.



Tube Extrusion Fixture measures continuous extrusion force of tube or sachet samples. Top vice grip included.



TA-MA (2)

Muco Adhesion Test Fixture Simulates body/temperature conditions and force needed to pull a tablet away from a mucosal surface.



TA-RT ②

Raft Tester measures alginate Raft Strength by pulling wire raft hook out of the sample material.



TA-MDI 🕑

Metered Dose Inhaler Fixture measures the push-button force to actuate the inhaler. Fixture Base Table required.



TA-STJ 📵

Syringe Test Fixture for measuring the force required to push or pull syringe plunger. Important to all syringe markets. Base Table required.



TA-TCA (2)

Tablet Coating Adhesion Fixture measures the adhesion force of a tablet coating to a tablet. Fixture Base Table required.



TA-BLS (2)

Bi-Layer Shear Fixture measures shear strength by cutting a two-part tablet or capsule using a guillotine blade. Fixture Base Table required.



TA-BPS (2)

Blister Pack Support Fixture is used to measure the force required to remove the tablet from its blister pack.
Fixture Base Table required.



TA-CLT (2)

Capsule Loop Tensile Test Fixture is used to measure the force required to split one half of a hard gel capsule. Fixture Base Table required.



TA-TCF ②

Tablet Compression Fixture performs a ring crush test (includes four rings). Bead compression test fixture is also included.

Fixture Base Table required.



TA-ATT 0 0

Adhesive Tack Tester for measuring stickiness of pressure sensitive adhesive materials such as tape. Rotary Base Table required.



TA-FSF 0 P

Film Support Fixture for puncture test to measure strength of fine films. Fixture Base Table required.



TA-LTT 🐠

Loop Tack Test measures the adhesive strength of pressure sensitive tape and stickers according to ASTM D6195.



TA-RIF (1)

Rotary Indexing Fixture measures adhesive force to pull tape off a surface. Multiple tests done on same sample for average value. Fixture Base Table required.



TA-GPJ 00

General Peeling Jig measures the adhesive strength needed to remove the lid from a sealed container at 0°, 45°, and 90° angles.



TA-PF180 (I)

180° Peel Fixture consists of top & bottom grips to measure adhesive strength when pulling tape off rigid surface using force at 180° angle.



TA-PF90 0

90° Peel Fixture measures the adhesive strength to pull a tape off of a rigid surface using force at a 90° angle.



TA-ATTPT **(I)** 

Adhesive Test Fixture measures adhesive force of plaster and similar samples. Includes flat cylindrical probe, base plate and top clamp to hold the sample.



TA-DGF (I) (D)

Dual Grip Fixture for tensile testing of thin films or integrity of seals for packaging.



TA-RCA 🐠 📵

Roller Cam Accessory grips measure the tensile strength and tear characteristics of material such as polymer films.



TA-SFJ 🐠 📵

Sliding Friction Jig measures the coefficient of friction for packaging materials according to ASTM D1894.



TA11 0

Cylindrical probe TA11 is used to measure the force vs. distance in mechanical springs. TA11: included in TA-P-KIT2 or as a standalone part.



### TA-P-KIT2 G @ D 00 D









General Probe Kit with Carrying Case

Variety of probes:

- (1) cutting wire
- (8) cylindrical
- (3) cone
- (1) knife-edge
- (2) ball
- (1) needle

Complies with BS and AACC.

# **Standard Probes**

TA5

Α	TA53	cutting wire
В	TA25/1000	50.8mm dia.
C	TA4/1000	38.1mm dia.
D	TA11/1000	25.4mm dia.
Ε	TA10	12.7mm dia.
F	TA17	24mm dia. 30°
G	TA15/1000	30mm dia. 45°
Н	TA2/1000	30mm dia. 60°
1	TA7	knife-edge
J	TA43	25.4mm dia.
K	TA18	12.7mm dia.
L	TA9	needle
M	TA41	6mm dia.
N	TA39	2mm dia.
0	TA44	4mm dia.

12.7mm dia.



# TA-P-KIT3 **6**

Curd Probe Kit one each of:

TA46 5mm TA47 8mm TA48 10mm



## **TA-CW-1500C**

Calibration Weight Set contains a combination of certified weights which may be used to confirm the calibration and linearity of each specific load cell.

### Also available:

A variety of cylindrical, cone, ball, needle, knife-edge, and wire cutting probes are also available.

AMETEK Brookfield can also make custom fixtures and probes for a variety of applications. Please contact AMETEK Brookfield or an authorized dealer to discuss solutions to your texture and materials testing challenges.

## **BASE TABLES**



#### TA-RT-KIT

Rotary Base Table provides quick and easy height adjustment to accommodate samples of various sizes. Included are pair of T-bolts for securing Rotary Base Table to slot in base of CT3 Texture Analyzer.



### **TA-BT-KIT**

Fixture Base Table is rectangular with removable insert which can be used as test surface or replaced with a number of test fixtures. Included are pair of T-bolts for securing Fixture Base Table to slot in base of CT3 Texture Analyzer. Also included are four sets of extension legs with different lengths to adjust the test surface height.



### **TA-ABT-KIT**

Adjustable Base Table for use with fixtures requiring table heights between 1.5" and 5"



The AMETEK Brookfield Texture Analysis Lab

# **TEXTURE ANALYSIS TESTING SERVICES**

A variety of texture analysis testing services can be performed at all AMETEK Brookfield locations (USA, UK, Germany, China and India). Most services are performed free-of-charge.

### **TEST AND RECOMMEND**

A simple evaluation designed to help determine the appropriate CT3 equipment for your application.

#### **SAMPLE PROFILING**

Analysis testing to determine specific properties of your sample.

### **DISPUTE RESOLUTION**

Analysis testing for mediating a resolution between producers and/or suppliers in cases where each has different results for the same material.

### **MULTIPLE SAMPLE TEST**

Expands the capability of your laboratory by utilizing AMETEK Brookfield's services to accomplish testing work on a timely basis.